

Produits Dentaires

Unique Quality
for Experts

Catalogue 8.4



GETZ HEALTHCARE DENTAL LIMITED

23/F, 8 Commercial Tower, 8 Sun Yip Street,
Chai Wan, Hong Kong SAR, People's Republic of China

香港柴灣新業街8號八號商業廣場23樓

Tel. : +852 2563 5193

Fax : +852 2565 9537

Email : info.dental.hk@getzhealthcare.com



Swiss quality
dental products



Swiss quality
dental products

Intro

Produits Dentaires SA is pleased to present the fourth edition of its new, updated, general catalogue N°8. This new edition is still divided into four main chapters, in order to facilitate the products search. The macro illustrations enable to enjoy the nature of each product as experienced by professionals in their daily practice.

Contact

Produits Dentaires SA
Rue des Bosquets 18
1800 Vevey
Switzerland
Phone +41 21 921 26 31
Fax +41 21 921 39 79
info@pdsa.ch
www.pdsa.ch



Produits Dentaires SA
Vevey . Switzerland

History

Founded in 1940 by William Gehrig, Produits Dentaires SA has advanced through the decades, thanks to reputed know-how, outstanding quality and innovative mindset. Today, the company is still family-owned and the third generation daily strives to implement actions answering the challenges of the international dental industry.

Social Responsibility

As a human-oriented and independent family business, Produits Dentaires SA has always considered staff integration a priority in its corporate culture. We also actively support local cultural events and international humanitarian associations through donations and free goods.

Products

Through this catalogue, Produits Dentaires SA offers to dentists a wide range of high quality medical devices used in endodontics, restoratives, prophylaxis, prosthetics and periodontology.

Ecological Commitment

Produits Dentaires SA is more than ever aware of the environmental issues related to the production and distribution of its products throughout the world. Compliance with requirements and constant optimisation of its processes show the willingness of the company to fulfill its duties towards our environment.

Distribution

Present in more than 100 countries on the 5 continents, Produits Dentaires SA has created, over the years, an extensive distribution network of agents, wholesale dealers and dental depots. Thanks to our partners, «PD» brand products have gradually gained international repute, synonymous with Swiss Quality at its best.

Innovation

Produits Dentaires SA invests in the development of new products which are based on the latest knowledge. As a Swiss manufacturer, we strive to present you always innovative novelties which meet our high quality standards and facilitate your daily work. In addition, existing products are continuously analyzed and improved. To meet this goal, Produits Dentaires SA works close together with a national and international expert team consisting of dentists, dental hygienists and other specialists from the medical field. In addition, several research projects are running in close cooperation with universities and colleges in Switzerland and worldwide.



Table of contents

| | | |
|----|-------------|------|
| 1. | Endodontics | 7-17 |
|----|-------------|------|

| | | |
|----|--------------|-------|
| 2. | Restoratives | 19-44 |
|----|--------------|-------|

| | | |
|----|---------------------------------|-------|
| 3. | Prophylaxis and Oral Hygiene | 45-48 |
|----|---------------------------------|-------|

| | | |
|----|--|-------|
| 4. | Special Supplies for the Clinic and the Laboratory | 51-56 |
|----|--|-------|

Available separately

ODUS Products
Medical Preparations

Orders Consult the Products Index at the end of the catalogue to find the relevant page where each item is clearly described and illustrated. Item code numbers are indicated under «REF» on a light gray background.



Endodontics

1.1 Points

| | |
|---------------------|---|
| Gutta-percha Points | 8 |
| Acrylic Points | 8 |
| Paper Points | 9 |

1.2 Instruments and Accessories

| | |
|------------------------------------|-------|
| Rotary Instruments | 10-11 |
| Map System | 13-15 |
| Endodontic Devices | 16 |
| Canal Irrigation and Desobturation | 17 |



Gutta-percha points

Precision points rolled and produced in compliance with the strictest hygienic requirements.
Radiopaque, Cadmium free, stiff yet flexible, of adequate plasticity for a successful condensation technique.

Regular type

Colour Pink
Shape Tapered
Length 28 mm



4 compartments box = 150 points

| Type | Ø | REF |
|---------|--------|-------|
| XX-fine | 0.5 mm | 23401 |
| X-fine | 0.6 mm | 23402 |
| Fine | 0.7 mm | 23403 |
| Medium | 1.0 mm | 23404 |

Standardized type

Head Flattened **Colour** Pink, colour coded
Sizes ISO # 15 to # 140 **Length** 29 mm

Box of 6 tubes = 300/150 points

| ISO # | | pieces | pink REF | colour coded REF |
|----------------|---|--------|-------------|------------------------|
| 15 | □ | 300 | 23412 | 23612 |
| 20 | ■ | 300 | 23413 | 23613 |
| 25 | ■ | 300 | 23414 | 23614 |
| 30 | ■ | 300 | 23415 | 23615 |
| 35 | ■ | 300 | 23416 | 23616 |
| 40 | ■ | 300 | 23417 | 23617 |
| 15-40 assort. | | 300 | 23418 | 23618 |
| 45 | □ | 300 | 23419 | |
| 50 | ■ | 300 | 23420 | |
| 55 | ■ | 300 | 23422 | |
| 60 | ■ | 300 | 23423 | |
| 70 | ■ | 300 | 23424 | |
| 80 | ■ | 300 | 23425 | |
| 45-80 assort. | | 300 | 23426 | |
| 90 | □ | 150 | 23430 | |
| 100 | ■ | 150 | 23431 | |
| 110 | ■ | 150 | 23432 | |
| 120 | ■ | 150 | 23434 | |
| 130 | ■ | 150 | 23435 | |
| 140 | ■ | 150 | 23436 | |
| 90-140 assort. | | 150 | 23437 | |



Acrylic points

Radiopaque and Cadmium free, they are mechanically produced to closest ISO tolerances and are stiff but flexible. They do not break or twist even in the narrowest and most curved canals and ensure a perfect obturation at the apical level. They finally ensure time-saving treatment as they do not need any condensation, though granting a tight apical obturation. They have the same rigidity as Silver points but can, as Gutta-percha points, be at any time eliminated from the root canal with the help of a bur or an organic solvent.

Disinfection

By immersion in a 2.5 % solution of Sodium Hypochlorite, rinsing with alcohol and drying.

Standardized type

Colours Pink, colour coded
Sizes ISO # 10 to # 140
Length 29 mm



Box of 6 tubes = 60 points

| ISO # | | pieces | pink REF | colour coded REF |
|----------------|------|--------|-------------|------------------------|
| 10 | 6×10 | 60 | 26702 | 26902 |
| 15 | 6×10 | 60 | 26703 | 26903 |
| 20 | 6×10 | 60 | 26704 | 26904 |
| 25 | 6×10 | 60 | 26705 | 26905 |
| 30 | 6×10 | 60 | 26706 | 26906 |
| 35 | 6×10 | 60 | 26707 | 26907 |
| 40 | 6×10 | 60 | 26708 | 26908 |
| 15-40 assort. | 6×10 | 60 | 26709 | 26909 |
| 45 | 6×10 | 60 | 26710 | 26910 |
| 50 | 6×10 | 60 | 26711 | 26911 |
| 55 | 6×10 | 60 | 26713 | 26913 |
| 60 | 6×10 | 60 | 26714 | 26914 |
| 70 | 6×10 | 60 | 26715 | 26915 |
| 80 | 6×10 | 60 | 26716 | 26916 |
| 45-80 assort. | 6×10 | 60 | 26717 | 26917 |
| 90 | 6×10 | 60 | 26718 | 26918 |
| 100 | 6×10 | 60 | 26719 | 26919 |
| 110 | 6×10 | 60 | 26720 | 26920 |
| 120 | 6×10 | 60 | 26722 | 26922 |
| 130 | 6×10 | 60 | 26723 | 26923 |
| 140 | 6×10 | 60 | 26724 | 26924 |
| 90-140 assort. | 6×10 | 60 | 26725 | 26925 |

Paper points

Rolled without binding agent, the «PD» Paper points are stiff but flexible. They correspond to the ISO requirements and are sterilized by gamma irradiation. Shelf life in the blister, 5 years.

Regular type

New Sterile blister
Colour White
Shape Tapered
Length 28 mm



6 compartments box = 200 points

| Type | Ø | REF |
|---------------------------|--------|-------|
| XX-fine | 0.7 mm | 25701 |
| X-fine | 0.8 mm | 25702 |
| Fine | 0.9 mm | 25703 |
| Medium | 1.1 mm | 25704 |
| Large | 1.3 mm | 25705 |
| X-large | 1.4 mm | 25706 |
| Assort. (XXF/XF/F/M/L/XL) | | 25707 |

Standardized type

New Sterile blister
Colours Plain
Sizes ISO # 15 to # 140
Length 29 mm



Box of 6 tubes = 300/240/120 points

| ISO # | | pieces | white REF |
|----------------|------|--------|-----------|
| 15 | 6×50 | 300 | 25780 |
| 20 | 6×50 | 300 | 25781 |
| 25 | 6×50 | 300 | 25782 |
| 30 | 6×50 | 300 | 25783 |
| 35 | 6×50 | 300 | 25784 |
| 40 | 6×50 | 300 | 25785 |
| 15-40 assort. | 6×50 | 300 | 25786 |
| 45 | 6×40 | 240 | 25787 |
| 50 | 6×40 | 240 | 25788 |
| 55 | 6×40 | 240 | 25789 |
| 60 | 6×40 | 240 | 25790 |
| 70 | 6×40 | 240 | 25791 |
| 80 | 6×40 | 240 | 25792 |
| 45-80 assort. | 6×40 | 240 | 25793 |
| 90 | 6×20 | 120 | 25794 |
| 100 | 6×20 | 120 | 25795 |
| 110 | 6×20 | 120 | 25796 |
| 120 | 6×20 | 120 | 25797 |
| 130 | 6×20 | 120 | 25798 |
| 140 | 6×20 | 120 | 25799 |
| 90-140 assort. | 6×20 | 120 | 25800 |

New Sterile blister **Colours** Colour coded
Sizes ISO # 15 to # 50 **Length** 29 mm



Box of 6 tubes = 300/240 points

| ISO # | | pieces | colour coded REF |
|---------------|--|--------|------------------|
| 15 | | 300 | 25860 |
| 20 | | 300 | 25861 |
| 25 | | 300 | 25862 |
| 30 | | 300 | 25863 |
| 35 | | 300 | 25864 |
| 40 | | 300 | 25865 |
| 15-40 assort. | | 300 | 25866 |
| 45 | | 240 | 25867 |
| 50 | | 240 | 25868 |

Gates-Glidden drills

Alloy Stainless Steel

Types RA, RA short

Lengths RA 32 mm, RA short 28 mm

Rotation speed 500-800 RPM



Rack of 6 pieces

| Type | N° | Ø | REF |
|----------|-----|------|-------|
| RA | 1 | 0.50 | 20302 |
| RA | 2 | 0.70 | 20312 |
| RA | 3 | 0.90 | 20322 |
| RA | 4 | 1.10 | 20332 |
| RA | 5 | 1.30 | 20342 |
| RA | 6 | 1.50 | 20352 |
| RA | 1-6 | – | 20362 |
| RA short | 1 | 0.50 | 20304 |
| RA short | 2 | 0.70 | 20314 |
| RA short | 3 | 0.90 | 20324 |
| RA short | 4 | 1.10 | 20334 |
| RA short | 5 | 1.30 | 20344 |
| RA short | 6 | 1.50 | 20354 |
| RA short | 1-6 | – | 20364 |

Peeso reamers

Alloy Stainless Steel

Types RA, RA short

Lengths RA 32 mm, RA short 28 mm

Rotation speed 800-1200 RPM



Rack of 6 pieces

| Type | N° | Ø | REF |
|----------|-----|------|-------|
| RA | 1 | 0.70 | 20402 |
| RA | 2 | 0.90 | 20412 |
| RA | 3 | 1.10 | 20422 |
| RA | 4 | 1.30 | 20432 |
| RA | 5 | 1.50 | 20442 |
| RA | 6 | 1.70 | 20452 |
| RA | 1-6 | – | 20462 |
| RA short | 1 | 0.70 | 20404 |
| RA short | 2 | 0.90 | 20414 |
| RA short | 3 | 1.10 | 20424 |
| RA short | 4 | 1.30 | 20434 |
| RA short | 5 | 1.50 | 20444 |
| RA short | 6 | 1.70 | 20454 |
| RA short | 1-6 | – | 20464 |

DRILL-A-POST . Spiral bladed tapered reamers

Alloy Stainless Steel

Type RA Shank

Point Blunt

Constant conicity 2°30'

Total length 31.5 mm

Working length 18.5 mm

Rotation speed 800 – 1200 RPM

Reamers corresponding to glass fiber posts FIBRAPOST PLUS and LUX N° 1 to 4 (page 40)

Packaging

Now also available in hygienic and resistant blisters of 4 reamers. Each reamer is packaged individually



Packaging classic Rack of 6 pieces

| N° | REF | N° | REF |
|-------------|-------|-------------|-------|
| 1 | 22501 | 1 | 22511 |
| 2 | 22502 | 2 | 22512 |
| 3 | 22503 | 3 | 22513 |
| 4 | 22504 | 4 | 22514 |
| 5 | 22505 | 5 | 22515 |
| 6 | 22506 | 6 | 22516 |
| 1-6 assort. | 22507 | 1-6 assort. | — |

Packaging blister Box of 4 pieces

| N° | REF |
|-------------|-------|
| 1 | 22651 |
| 2 | 22652 |
| 3 | 22653 |
| 4 | 22654 |
| 1-4 assort. | 22660 |

Root canal barbed broaches

Alloy Stainless Steel blade

Length 21 mm

Spiral handle Colour coded



Rack of 12 pieces

| ISO # | colour code | REF |
|---------------|-------------|-------|
| 25 | xxx-fine | 21151 |
| 30 | xx-fine | 21152 |
| 35 | x-fine | 21153 |
| 40 | fine | 21154 |
| 50 | medium | 21155 |
| 60 | coarse | 21156 |
| 25-60 assort. | | 21157 |





MAP SYSTEM

The MAP (Micro-Apical Placement) System provides a unique and efficient method for placing root canal repair materials.

- either by orthograde obturation for the treatment of perforations, root-end fillings and pulp-cappings using curved, straight or NiTi « memory shape » needles.
- or by retrograde obturation after apical resection thanks to the specially designed, triple-angled needles (left and right angled) and/or hooked needles.

1.2 MAP SYSTEM

| | |
|-------------------------|----|
| Kits | 14 |
| PD MTA White | 15 |
| Needles and Accessories | 15 |

Kits



Intro kit NiTi

REF 20294

For orthograde and retrograde obturations of anterior teeth

Contains

| | |
|---------------------------------------|------------------|
| 1 sterilisable aluminium case | 185 × 95 × 30 mm |
| 1 stainless syringe | |
| 1 NiTi Memory Shape needle | N° 0 (yellow) |
| + 3 spare PEEK plungers | N° 0 |
| + 1 cleaning curette | N° 0 |
| 1 curved needle (incl. NiTi plunger) | N° 0 |
| 1 NiTi Memory Shape needle | N° 1 (red) |
| + 3 spare PEEK plungers | N° 1 |
| + 1 cleaning curette | N° 1 |
| 1 curved needle (incl. NiTi plunger) | N° 1 |
| 1 stainless steel well | ø 14 mm |
| 1 paste cleaning device (key + stops) | N° 0 |
| 1 paste cleaning device (key + stops) | N° 1 |



Universal kit NiTi

REF 20295

For orthograde and retrograde obturations of anterior and posterior teeth

Contains

| | |
|---------------------------------------|-------------------|
| 1 sterilisable aluminium case | 185 × 145 × 30 mm |
| 1 stainless syringe | |
| 1 NiTi Memory Shape needle | N° 0 (yellow) |
| 1 Hooked needle | N° 0 |
| + 6 spare PEEK plungers | N° 0 |
| + 1 cleaning curette | N° 0 |
| 1 NiTi Memory Shape needle | N° 1 (red) |
| 1 Hooked needle | N° 1 |
| + 6 spare PEEK plungers | N° 1 |
| + 1 cleaning curette | N° 1 |
| 1 curved needle (incl. NiTi plunger) | N° 0 |
| 1 curved needle (incl. NiTi plunger) | N° 1 |
| 1 stainless steel well | ø 14 mm |
| 1 paste cleaning device (key + stops) | N° 0 |
| 1 paste cleaning device (key + stops) | N° 1 |



Surgical kit

REF 20270

For retrograde obturations

Contains

| | |
|---------------------------------------|-------------------|
| 1 sterilisable aluminium case | 185 × 145 × 30 mm |
| 1 stainless syringe | |
| 1 triple-angle needle, right | N° 0/R (yellow) |
| 1 triple-angle needle, left | N° 0/L |
| 1 triple-angle needle, «retro» right | N° 0/R |
| 1 triple-angle needle, «retro» left | N° 0/L |
| + 6 spare PEEK plungers | N° 0 |
| + 1 cleaning curette | N° 0 |
| 1 triple-angle needle, right | N° 1/R (red) |
| 1 triple-angle needle, left | N° 1/L |
| + 6 pistons PEEK | N° 1 |
| + 1 cleaning curette | N° 1 |
| 1 stainless steel well | ø 14 mm |
| 1 paste cleaning device (key + stops) | N° 0 |
| 1 paste cleaning device (key + stops) | N° 1 |

PD MTA White

Endodontic filling material, specially developed to be placed with MAP System

Advantages

- Optimized particle size
- Avoids bacterial migration
- Excellent marginal sealing capacity
- Stimulates the formation of a dentine layer (pulp capping)
- 2 sachets with 280 mg each



Packaging

2 applications (2 x 280 mg) and 1 bottle of 3 ml pure H₂O

REF
20190

Needles



NiTi «memory shape» needle
(2 PEEK plungers included)

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20275 |
| 1 | red | 1.10 | 20277 |
| 2 | blue | 1.30 | 20278 |



Straight needle (1 NiTi plunger included)

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20238 |
| 1 | red | 1.10 | 20239 |
| 2 | blue | 1.30 | 20292 |
| 3 | green | 1.50 | 20293 |



Curved needle (1 NiTi plunger included)

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20236 |
| 1 | red | 1.10 | 20237 |
| 2 | blue | 1.30 | 20290 |
| 3 | green | 1.50 | 20291 |



Triple-angle needle (1 PEEK plunger included)

| N° | type | colour | ext. Ø mm | REF |
|----|-------|--------|-----------|-------|
| 0 | right | yellow | 0.90 | 20232 |
| 0 | left | yellow | 0.90 | 20233 |
| 1 | right | red | 1.10 | 20234 |
| 1 | left | red | 1.10 | 20235 |



Triple-angle «retro» needle
(1 PEEK plunger included)

| N° | type | colour | ext. Ø mm | REF |
|----|-------|--------|-----------|-------|
| 0 | right | yellow | 0.90 | 20271 |
| 0 | left | yellow | 0.90 | 20272 |
| 1 | right | red | 1.10 | 20273 |
| 1 | left | red | 1.10 | 20274 |



Hooked needle (1 PEEK plunger included)

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20263 |
| 1 | red | 1.10 | 20264 |
| 2 | blue | 1.30 | 20265 |
| 3 | green | 1.50 | 20266 |

Accessories



Plunger NiTi

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20247 |
| 1 | red | 1.10 | 20248 |
| 2 | blue | 1.30 | 20283 |
| 3 | green | 1.50 | 20284 |



Plunger PEEK (6 pieces)

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20245 |
| 1 | red | 1.10 | 20246 |
| 2 | blue | 1.30 | 20281 |
| 3 | green | 1.50 | 20282 |



Cleaning curette NiTi

| N° | colour | ext. Ø mm | REF |
|----|--------|-----------|-------|
| 0 | yellow | 0.90 | 20240 |
| 1 | red | 1.10 | 20241 |
| 2 | blue | 1.30 | 20243 |
| 3 | green | 1.50 | 20244 |



Paste cleaning device

| N° | colour | REF |
|----|--------|-------|
| 0 | yellow | 20252 |
| 1 | red | 20253 |
| 2 | blue | 20254 |
| 3 | green | 20255 |

« Endo » plastic pipettes

For root canal or cavity filling with fluids.

Capacity ~0.3 ml

Points PE Fine, flexible and unbreakable PE tip, acid / solvent resistant tips

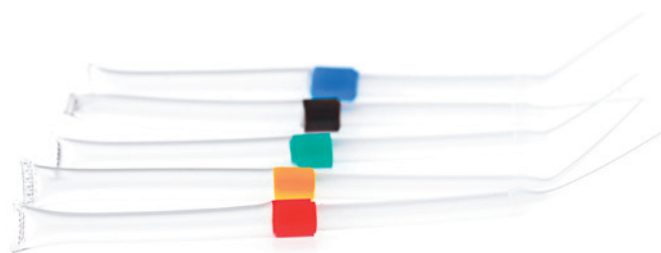
Transparent

Colour coded Red, blue, yellow, green and black

Sterilizable Ethylene Oxide gas

(do not use heat autoclave)

Desinfection The pipettes may also be disinfected in a 3% solution of Sodium Hypochlorite, in Hydrogen Peroxide or in Alcohol



Box

5 pipettes of assorted colours

REF

11750

« One-piece » plastic pipettes

For inserting fluids into the narrowest canals.

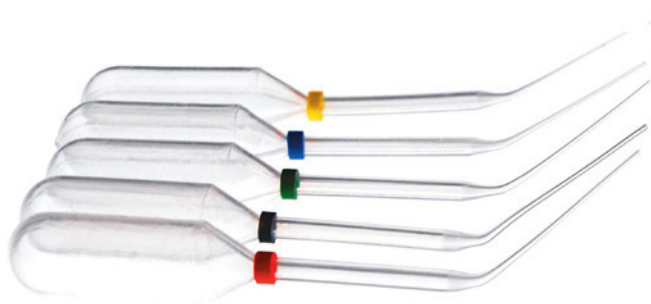
Capacity ~1 ml

Points PE Ultra-fine tips, flexible and resistant, acid and solvent resistant

Transparent

Colour coded Red, blue, yellow, green and black

Desinfection The pipettes may also be disinfected in a 3% solution of Sodium Hypochlorite, in Hydrogen Peroxide or in Alcohol



Box

5 pipettes of assorted colours

REF

11752

Marking rings

For the identification of hand instruments.



Box of 48 pieces

| Colour | Ø 6,8 mm REF | Ø 8 mm REF | Ø 10 mm REF |
|---------|-----------------|---------------|----------------|
| □ | 21341 | 21361 | 21381 |
| ■ | 21342 | 21362 | 21382 |
| ■ | 21343 | 21363 | 21383 |
| ■ | 21344 | 21364 | 21384 |
| ■ | 21345 | 21365 | 21385 |
| ■ | 21346 | 21366 | 21386 |
| ■ | 21347 | 21367 | 21387 |
| ■ | 21348 | 21368 | 21388 |
| ■ | 21349 | 21369 | 21389 |
| ■ | 21350 | 21370 | 21390 |
| Assort. | 21351 | 21371 | 21391 |

EDTA 17 %

Solution for root canal irrigation.
Gel for the chemical widening of the root canal.

Properties . Indication

EDTA 17 % is a solution/gel to chemically facilitate the preparation of the root canal during its instrumental shaping (gel) and irrigation (solution). It acts by chelation of the Calcium ions and elimination of inorganic debris (smear layers) through the softening of the dentinal walls.

Composition

1 ml of solution / gel contains:
Edetate tetrasodic (EDTA) 190 mg (17% w/w),
cetyltrimethylammonium Bromide 10 mg,
excipient ad 1 ml



Packaging

| | REF |
|----------------------------|-------|
| Bottle of 15 ml (solution) | 11565 |
| Syringe of 3 g (gel) | 11566 |

Eugenate Desobturator

Solution for removing oxy-eugenolates and Gutta-percha obturations from the root canal.

Properties . Indication

Eugenate Desobturator is a solution to be used specifically for the softening and subsequent mechanical removal of eugenate type cements and Gutta-percha obturations from the root canal.

Composition

Oil of Orange, excipient



Packaging

| | REF |
|-----------------|-------|
| Bottle of 15 ml | 11555 |

Resin Remover

Solution for removing resin-based fillings from the root canal.

Properties . Indication

Resin Remover is a solution to be used specifically for the softening and subsequent mechanical removal of phenoplastic resins from the root canal.


Composition

N,N-Dimethylformamide 15 %,
Oil of Citronella 15 %, Tincture
of Benzoin 5 %, Tolu Balsam 5 %,
Essential Oils 4.6 %, excipient ad 100 %



Packaging

| | REF |
|-----------------|-------|
| Bottle of 15 ml | 11556 |

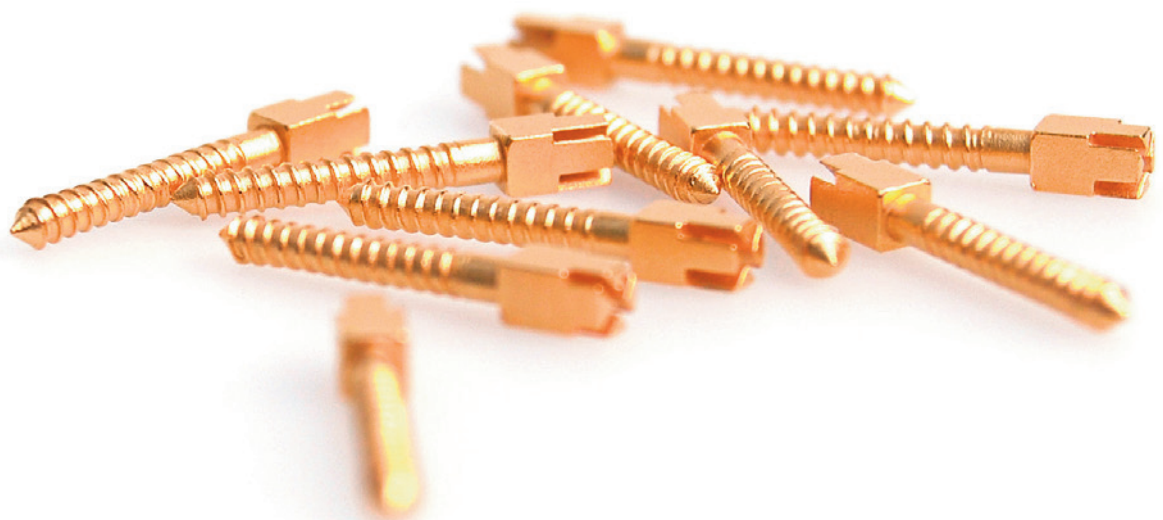

Restoratives

2.1 Materials and Preparations

| | |
|-----------------------------|-------|
| Liners and Varnishes | 20-22 |
| Temporary Filling Materials | 23 |
| Cements | 24-27 |

2.2 Devices and Instruments

| | |
|---|-------|
| Crown Forms, Rings and Crowns | 28-32 |
| Matrices, Wedges, Articulating Paper and Adhesive | 33-35 |
| Abrasives | 36-37 |
| FIBRAPOST & SEALACORE | 39-41 |
| Screw Posts | 42-43 |



«Cavity Liner» paste

Calcium hydroxide, ready-to-use, radiopaque, pre-mixed paste.

Properties

Calcium Hydroxide Paste «Cavity Liner» is a creamy paste used as a liner or as a temporary dressing. Due to its alkaline pH (~12), Calcium Hydroxide neutralizes the free acids of the obturating cements and of the restorative materials. It reduces the hypersensitivity of the dentine and helps the neoformation of secondary dentine.

Indications

As a liner

- ▶ Pulp capping on vital pulp (indirect capping)
- ▶ Capping after pulpotomy (direct capping)
- ▶ Cavity lining before filling

As a temporary dressing

- ▶ Root canal treatment and periapical therapy

Composition

1 ml of paste contains 540 mg Calcium Hydroxide and 70 mg Barium Sulfate in a Methylcellulose base



Packaging

Syringe of 2.5 ml + 2 applicator tips

REF

11020

«Cavity Liner» light cure

Calcium hydroxide, light-curing, radiopaque paste for cavity lining.

Properties

Radiopaque light-curing, one-component paste based on Calcium Hydroxide and Dimethacrylates

This paste provides:

- ▶ Excellent stability
- ▶ High resistance to resorption and dissolution
- ▶ Excellent adhesion to various materials

Indications

- ▶ Cavity lining under composites or amalgams
- ▶ Indirect pulp capping

Composition

Urethane Dimethacrylate, Calcium Hydroxide, Barium Sulfate, Silicates, excipient



Packaging

Syringe of 4 g + 5 nozzles

REF

11035

« Cavity Liner » compound

Two-components (self-curing) radiopaque Calcium Hydroxide composition.

Properties

The Calcium Hydroxide compound is a two component paste, of ivory colour, used as a liner or a pulp capper under all usual filling materials

Indications

- Pulp capping on vital pulp (indirect capping)
- Capping after pulpotomy (direct capping)
- Cavity lining before filling

Precautions

- Due to its alkalinity the Calcium Hydroxide Compound should only be used where it can be protected from the intra-oral environment
- It should not be used as a temporary filling

Composition

Catalyst: Calcium Hydroxide, Polymeric plasticizer, Zinc Oxide

Base: Salicylic Acid derivatives, Calcium Wolframate, Calcium Sulfate



Packaging

1 tube base 13 g
+ 1 tube catalyst 11 g
+ mixing pads

REF

11030

Calcium Hydroxide powder

Chemically pure Calcium Hydroxide in powder form.



Box

50 g

REF

11025

Copal varnish

Protective varnish for the insulation of exposed dentine.

Properties . Indications

Copal varnish is the ideal cavity liner and dentinal tubuli sealer. It ensures perfect protection to the exposed surfaces of dentine in teeth prepared for permanent fillings or restorations and is useful as an intermediary varnish prior to the insertion of amalgams, cements or composites.

Composition

Varnish Gum Damar 2.4 %, Resin, Linseed Oil derivatives, Chlorobutanol 0.1 %, Solvent organic ad 100 %

Solvent Ethyl Acetate

Note These products are subject to transport restrictions



Packaging

| | | REF |
|-------------------|-------------------|-------|
| Varnish + solvent | bottle 15 ml each | 11040 |
| Varnish alone | bottle 15 ml | 11041 |

Vaseline

Pure white Petrolatum for dental use.



Tube

| | REF |
|-------|-------|
| ~15 g | 10651 |

TEMPOTEC . Temporary Eugenol free filling material in paste

Ready to use, self-setting filling material for temporary fillings.

Properties

TEMPOTEC is a ready to use, slightly radiopaque paste. Its excellent adhesion to the dentine and its sealing qualities ensure a hermetic filling. The cement can be applied to the cavity whether it is moist or dry and causes no irritation to either the pulp or the gums. It hardens after a few minutes in presence of moisture (saliva) and can easily be removed with a probe.

Indications

- ▶ Temporary filling of cavities
- ▶ Temporary sealing on top of medicaments
- ▶ Temporary sealing of posts

Composition

Zinc Oxide , Zinc sulfate, Calcium sulfate, excipient



Packaging

| | | REF |
|-----|------|-------|
| Jar | 40 g | 10204 |

Gutta-percha stopping

Radiopaque, temporary, Gutta-percha based, filling material.

Sticks of about ~Ø 4 mm, length 100 mm

Melting point ~45°C

Colours White, pink, ivory or assorted

Composition Gutta-percha 11.3%, Calcium Carbonate 9.6%, Zinc Oxide 59.7%, excipient ad 100%



| Colour | 30 g REF | 120 g REF |
|----------|-------------|--------------|
| White | 10223 | 10228 |
| Pink | 10222 | 10227 |
| Ivory | 10221 | 10226 |
| Assorted | 10220 | 10225 |

Pro-Pul-Pan . Zinc Oxide-Eugenol reinforced cement

Temporary, Polymer reinforced root canal cement.

Properties

The therapeutic value of Zinc Oxide-Eugenol cement Pro-Pul-Pan is based on its compatibility with both hard and soft tissues of the mouth as well as on its slight topical effect. Zinc Oxide-Eugenol cements, like Pro-Pul-Pan, are particularly reputed for their sedative effect upon pulpal tissues, especially when restoring teeth with deep carious lesions. Additionally, it provides a good marginal seal and thermal insulation

Indications

- Temporary fillings, thermal insulation, pulp capping
- Periodontal dressings
- Root canal sealer
- Temporary sealing of crowns and inlays

Colour White

Composition

Powder Zinc Oxide, Zinc Acetate, Polymethylmetacrylate
Liquid Eugenol



| Packaging | | REF |
|-----------|----------------------------|-------|
| Portion | 30 g powder + 15 ml liquid | 10260 |
| Powder | bottle 30 g | 10261 |

C&B Zinc Phosphate cement

Luting cement for crowns, bridges and inlays.

Properties . Indications

A universal, all-purpose cement for cementing crowns, bridges, inlays, orthodontic appliances, etc. It can also be used as a filling material, a cavity liner and for pulp capping under gold, amalgam, porcelain and acrylic fillings

Shades

- N° 1 cream white
- N° 2 yellow white

Composition

Powder Zinc Oxide, Magnesium Oxide, Aluminium Phosphate

Liquid Phosphoric Acid, Aluminium Hydroxide, excipient



| Packaging | N° 1 REF | N° 2 REF |
|-----------|------------------------------------|-------------|
| | Portion 30 g powder + 15 ml liquid | 10230 |
| | | 10233 |

Rosen root canal cement . Rickert's paste

Radiopaque root canal sealer.

Properties

Mixed to a paste, this cement provides a very suitable material for permanent root canal sealing and filling. It fulfills practically every criteria as an ideal sealer. It is non-irritant, dimensionally stable, radiopaque and does not dissolve in the tissue fluids. Moreover it adapts well to the walls of the canal, is easy to shape and to remove.

Caution

Due to the presence of precipitated Silver, this sealer can stain and should only be used in posterior teeth

Indications

- Permanent sealing of the root canal
- Endodontic sealer which can be used with all types of endodontic points (Gutta-percha, Acrylic, Silver or Titanium)

Composition

Powder Zinc Oxide 41.2%, Thymol Iodide 12.8%, Rosin 16%, Precipitated Silver 30 %

Liquid Oil of Clove 78 %, Canada Balsam 22 %



Packaging

Portion 15 g powder + 15 ml liquid

REF

10250

Rosen root canal cement Silver free

Radiopaque Root canal sealer without Silver.

Properties

This sealer is manufactured according to Rickert's formula (see above) but does not include precipitated Silver. This avoids any risk of staining the teeth but, on the other hand, also reduces the bactericidal action of the paste.

Indications

- Permanent sealing of the root canal
- Sealer which can be used with all types of endodontic points (Gutta-percha, Acrylic, Silver or Titanium)

Composition

Powder Zinc Oxide 41.2%, Thymol Iodide 12.8%, Rosin 16%, Radiopaque excipient ad 100 %

Liquid Oil of Cloves 78 %, Canada Balsam 22 %



Packaging

Portion 15 g powder + 15 ml liquid

REF

10245

Glass Ionomer cements



Ionover Filling

Glass Ionomer Filling cement

Properties

Glass Ionomer Filling cement is radiopaque and offers excellent adhesion to both enamel and dentine. As it prevents the formation of marginal leakage, retentions are not strictly necessary though undercuts provide additional hold to the filling. The thickness of the filling layer should be at least 1 mm

Indications

- Restoration (especially class I) of deciduous teeth
- Repair of conical lesions and enamel defects as well as root erosions in the dental neck area, including class V fillings
- Obturation of cavities with carious lesions
- Class III fillings and small class I fillings, including extended fissure sealing
- Core build-up under crowns

Shades

C1 universal
B1 light

Composition

Powder Fluoro-Phospho-Silicate of Calcium and Aluminium
Liquid Polyacrylic acid

Packaging

Portion 15 g powder with a measuring scoop
+10 ml of liquid with a drop dispenser

C1
REF

10390

B1
REF

10392



Ionover Luting

Glass Ionomer Luting cement

Properties

Glass Ionomer Sealing (Luting) cement offers excellent adhesion to both enamel and dentine as well as to metallic surfaces. Its low level of solubility in the mouth and low acidity are its chief assets

Indications

- Sealing of crowns, bridges, inlays and onlays
- Sealing of crown posts
- Sealing of orthodontic rings

Shade

White

Composition

Powder Fluoro-Phospho-Silicate of Calcium and Aluminium
Liquid Polyacrylic acid

Packaging

Portion 15 g powder with a measuring scoop
+10 ml of liquid with a drop dispenser

REF

10385

Zinc Oxide

Mixed with Eugenol to a paste, this powder becomes a temporary cement.

Indications Temporary obturation of the root canal Protection of exposed pulp

Composition Zinc Oxide > 99.9%



| Jar | REF |
|-------|-------|
| 100g | 11431 |
| 250g | 11432 |
| 500g | 11430 |
| 1000g | 11433 |

Eugenol

Essential oil, obtained from purified oil of Cloves.

Properties

Mixed with Zinc Oxide to a paste, it becomes a temporary cement. It is also used for pulp capping.

Indications

- Temporary obturation of the root canal
- Protection of exposed pulp

Composition

Pure Eugenol (>99%)



| Bottle | REF |
|--------|-------|
| 15ml | 11420 |
| 30ml | 11421 |
| 100ml | 11422 |

Jiffy tubes for cement

Properties

Transparent, resistant and flexible acetate of cellulose disposable tubes for the safe and easy introduction of medicaments and/or materials of creamy consistency into cavities or into the pulp chamber. They are also widely used in the crown form technique where PELLA Transparent crown forms are used as a matrix (page 28)


Shapes Straight or curved



| Packaging | REF |
|------------------------------------|-------|
| Straight box of 36 tubes | 10681 |
| Curved box of 24 tubes | 10680 |
| Assortment 18 straight + 12 curved | 10682 |

PELLA transparent crown forms

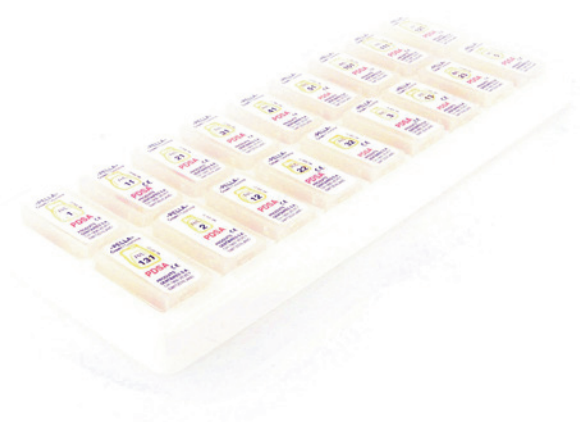- For the «peel-off» technique
- 78 different shapes, namely 39 pairs right (R) and left (L), reproducing the original form of the tooth
- Right crowns are identified in yellow colour and left crowns in pink













| Plexi-box | REF |
|----------------------|-------|
| 78 crowns (39 pairs) | 33410 |

| Assortment of 18 pairs | 33405 |
|---|-------|
| RL 1, 11, 21, 31, 41, 51, 101, 111, 121, 131, 2, 12, 22, 32, 3, 13, 23, 5 | |



Refill: 6-pack box = 6 pairs

| R/L | REF | R/L | REF |
|-----|-------|-----|-------|
| 1 | 33421 | 33 | 33621 |
| 2 | 33431 | 41 | 33631 |
| 3 | 33441 | 42 | 33641 |
| 4 | 33451 | 43 | 33651 |
| 5 | 33461 | 51 | 33661 |
| 6 | 33471 | 101 | 33671 |
| 7 | 33481 | 111 | 33681 |
| 11 | 33491 | 121 | 33691 |
| 12 | 33501 | 131 | 33701 |
| 13 | 33511 | 141 | 33711 |
| 14 | 33521 | 151 | 33721 |
| 15 | 33531 | 161 | 33731 |
| 16 | 33541 | 171 | 33741 |
| 17 | 33551 | 201 | 33751 |
| 21 | 33561 | 202 | 33761 |
| 22 | 33571 | 203 | 33771 |
| 23 | 33581 | 211 | 33781 |
| 24 | 33591 | 212 | 33791 |
| 31 | 33601 | 213 | 33811 |
| 32 | 33611 | | |

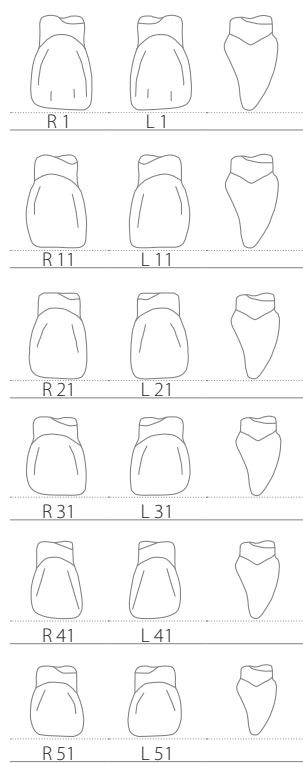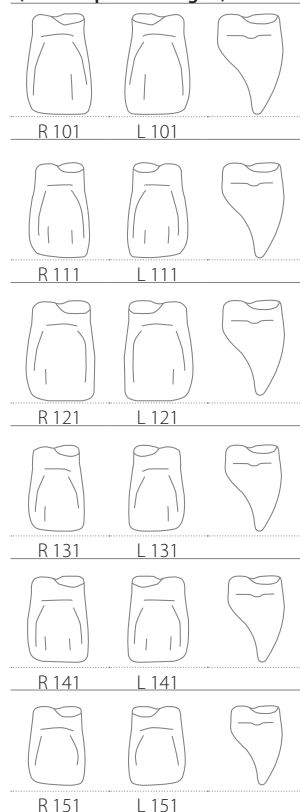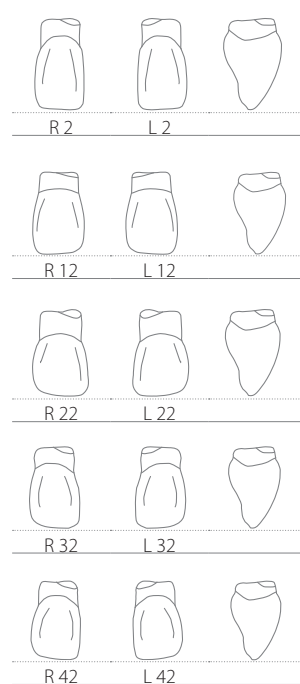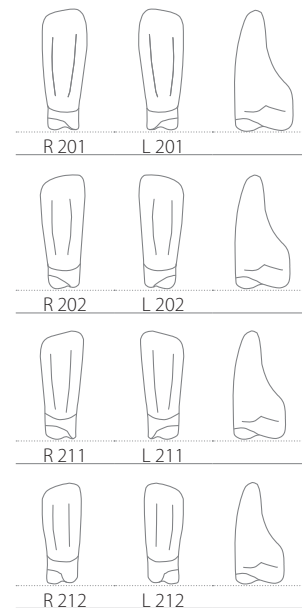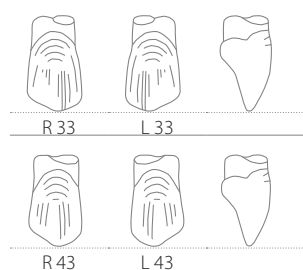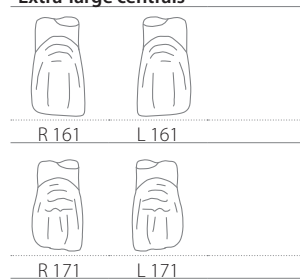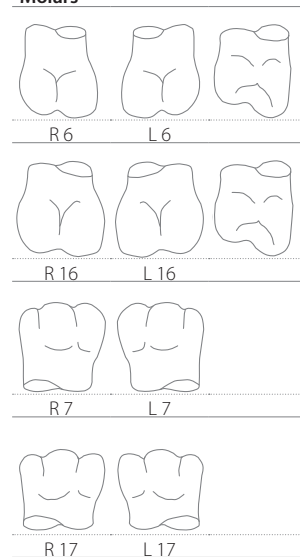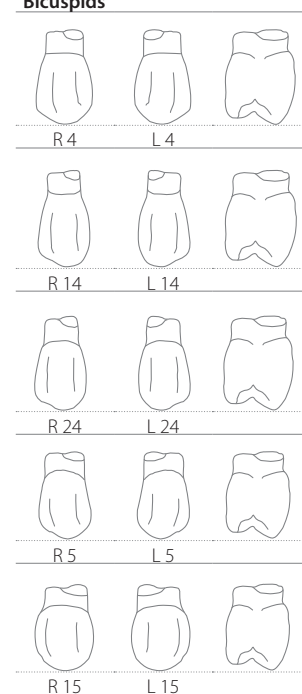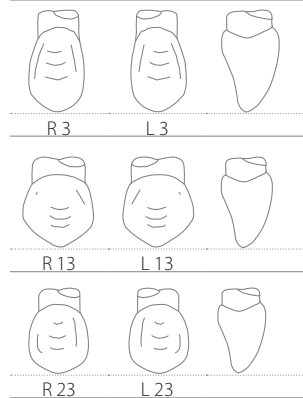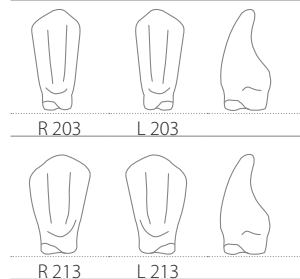**Upper Centrals
(normal forms)****Upper Centrals
(cervical parts enlarged)****Upper
laterals****Lower
incisors****Upper Extra-large
canines****Extra-large centrals****Molars****Bicuspids****Upper canines****Lower canines**

Copper bands

For impression of inlays by the indirect technique.

Types Hard or soft

24 differents Stamped size (circumference)

Height 12.5 mm

Thickness Wall 0.15 mm



Box assortment of 100 pieces

N° 1-20 (5 pieces each size, ½ sizes no included)

| | |
|------|-------|
| Hard | 34100 |
| Soft | 34101 |

Refill of 12 pieces

| N° | circumference (mm) int. | int. Ø | hard REF | soft REF |
|-----|----------------------------|--------|-------------|-------------|
| 1 | 15 | 4.8 | 34112 | 34110 |
| 2 | 16 | 5.1 | 34122 | 34120 |
| 3 | 17 | 5.4 | 34132 | 34130 |
| 4 | 18 | 5.7 | 34142 | 34140 |
| 4½ | 19 | 6.1 | 34152 | 34150 |
| 5 | 20 | 6.4 | 34162 | 34160 |
| 6 | 21 | 6.7 | 34172 | 34170 |
| 7 | 22 | 7.0 | 34182 | 34180 |
| 8 | 23 | 7.3 | 34192 | 34190 |
| 8½ | 24 | 7.7 | 34202 | 34200 |
| 9 | 25 | 8.0 | 34212 | 34210 |
| 10 | 26 | 8.3 | 34222 | 34220 |
| 11 | 27 | 8.6 | 34232 | 34230 |
| 12 | 28 | 8.9 | 34242 | 34240 |
| 12½ | 29 | 9.2 | 34252 | 34250 |
| 13 | 30 | 9.6 | 34262 | 34260 |
| 14 | 31 | 9.9 | 34272 | 34270 |
| 15 | 32 | 10.2 | 34282 | 34280 |
| 16 | 33 | 10.5 | 34292 | 34290 |
| 16½ | 34 | 10.8 | 34302 | 34300 |
| 17 | 35 | 11.1 | 34312 | 34310 |
| 18 | 36 | 11.5 | 34322 | 34320 |
| 19 | 37 | 11.8 | 34332 | 34330 |
| 20 | 38 | 12.1 | 34342 | 34340 |

Transparent rings

Properties

Transparent, flexible, thermoplastic, and resistant acetate cellulose rings. After softening in water heated at about 50°C, their thermoplastic quality ensures adaptation to the cervical outline.

Height 9 mm

Transparent rings can be used :

1. for core-build-up
2. as matrices with all kinds of dental filling, particularly with light-curing materials, like composites



Box assortment of 100 pieces (5 pieces each size)
1-20 **10760**

Refill of 12 pieces

| N° | circum. int. | REF |
|----|--------------|-------|
| 1 | 15 | 10701 |
| 2 | 16 | 10702 |
| 3 | 17 | 10703 |
| 4 | 18 | 10704 |
| 5 | 20 | 10706 |
| 6 | 21 | 10707 |
| 7 | 22 | 10708 |
| 8 | 23 | 10709 |
| 9 | 25 | 10711 |
| 10 | 26 | 10712 |
| 11 | 27 | 10713 |
| 12 | 28 | 10714 |
| 13 | 30 | 10716 |
| 14 | 31 | 10717 |
| 15 | 32 | 10718 |
| 16 | 33 | 10719 |
| 17 | 35 | 10721 |
| 18 | 36 | 10722 |
| 19 | 37 | 10723 |
| 20 | 38 | 10724 |

Preformed Aluminium crowns

Prosthetic temporary crown during endo-treatment.

Types Hard or soft
20 sizes Stamped
Wall Thickness 0.3 mm



Box assortment of 100 pieces N° 1-20 (5 pieces each size)

| | |
|------|-------|
| Hard | 36100 |
| Soft | 36101 |

Bicuspid upper refill of 12 pieces

| N° | A | B | C | hard REF | soft REF |
|----|-----|-----|------|-------------|-------------|
| 1 | 7.4 | 6.5 | 8.2 | 36112 | 36110 |
| 2 | 7.6 | 6.8 | 8.6 | 36122 | 36120 |
| 3 | 7.6 | 7.8 | 9.2 | 36132 | 36130 |
| 4 | 8.0 | 8.2 | 9.6 | 36142 | 36140 |
| 5 | 7.6 | 8.6 | 10.0 | 36162 | 36160 |

Bicuspid lower refill of 12 pieces

| | | | | | |
|----|-----|-----|-----|-------|-------|
| 6 | 8.0 | 6.2 | 7.3 | 36172 | 36170 |
| 7 | 8.2 | 6.5 | 7.7 | 36182 | 36180 |
| 8 | 8.2 | 6.9 | 8.0 | 36192 | 36190 |
| 9 | 7.9 | 7.2 | 8.3 | 36212 | 36210 |
| 10 | 7.6 | 7.9 | 8.6 | 36222 | 36220 |


Molars lower refill of 12 pieces

| | | | | | |
|----|-----|------|------|-------|-------|
| 11 | 7.0 | 9.8 | 9.2 | 36232 | 36230 |
| 12 | 7.9 | 10.3 | 9.6 | 36242 | 36240 |
| 13 | 7.7 | 10.8 | 9.9 | 36262 | 36260 |
| 14 | 7.9 | 11.2 | 10.2 | 36272 | 36270 |
| 15 | 9.5 | 11.5 | 10.5 | 36282 | 36280 |

Molars upper refill of 12 pieces

| | | | | | |
|----|-----|------|------|-------|-------|
| 16 | 8.2 | 9.6 | 10.4 | 36292 | 36290 |
| 17 | 8.9 | 10.8 | 11.1 | 36312 | 36310 |
| 18 | 8.5 | 10.9 | 11.4 | 36322 | 36320 |
| 19 | 9.1 | 11.3 | 11.9 | 36332 | 36330 |
| 20 | 8.8 | 11.6 | 12.0 | 36342 | 36340 |

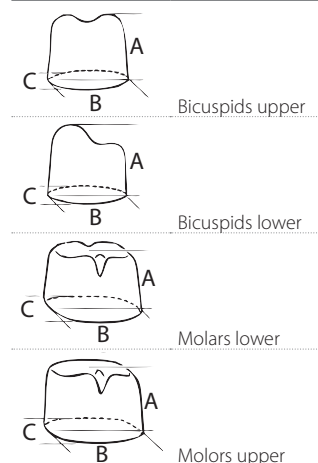Sizes



Tofflemire matrix retainers

Stainless Steel, designed as per Dr Tofflemire's original model.



Box with matrix retainer + 12 bands


| Shape | Content | REF |
|-----------|--------------------------|-------|
| Junior | 12 assorted bands . thin | 41501 |
| Senior | 12 assorted bands . thin | 41502 |
| Universal | 12 assorted bands . thin | 41503 |
| Junior | retainer alone | 41504 |
| Senior | retainer alone | 41505 |
| Universal | retainer alone | 41506 |

Tofflemire matrix bands

Alloy Stainless Steel

Thickness Thin 0.05 mm, x-thin 0.035 mm

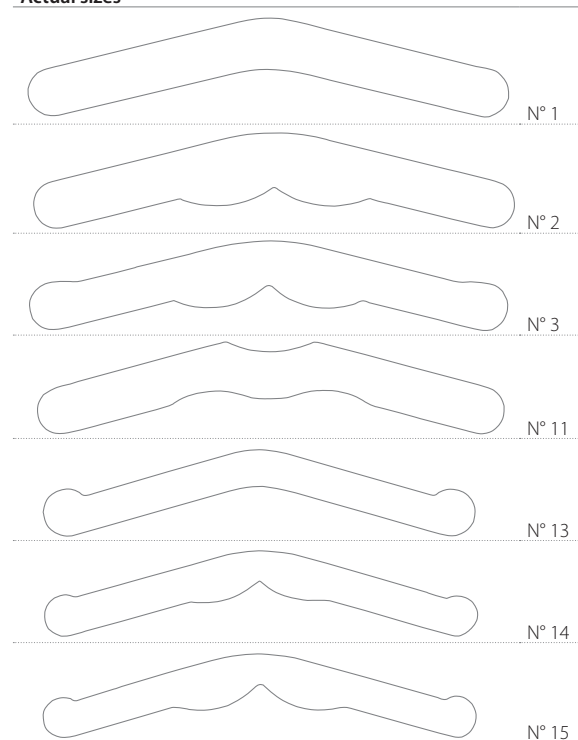Actual sizes



Pack of 12 pieces

| Adults | thin REF | x-thin REF |
|--------------------------|-------------|---------------|
| # 1 | 41521 | 41511 |
| # 2 | 41522 | 41512 |
| # 3 | 41523 | 41513 |
| # 11 | 42524 | 41514 |
| Children | | |
| # 13 | 41525 | 41515 |
| # 14 | 41526 | 41516 |
| # x15 | 41527 | 41517 |
| Assorted | | |
| 4× #1 – 3× #2 – 1× #3-15 | 41528 | 41518 |

Interdental Wooden wedges





Wood Birchwood

Types Plain, colour coded

Sizes Small, medium, long, wide



Box of 100 and 400 pieces

| Size | 100 pcs REF | 400 pcs REF | colour coded | 100 pcs REF | 400 pcs REF |
|----------|----------------|----------------|--|----------------|----------------|
| Small | 21721 | 21731 | plain  | 21741 | — |
| Medium | 21722 | 21732 | red  | 21742 | — |
| Long | 21723 | 21733 | yellow  | 21743 | — |
| Wide | 21724 | 21734 | blue  | 21744 | — |
| Assorted | 21725 | 21735 | assorted | 21745 | 21755 |

Anatomical Wooden wedges

Anatomically shaped Interdental

Maplewood wedges for a perfect point of contact





Square extremity for a tight grip

Pointed edge for safe, atraumatic insertion

Type Colour coded



Box of 100 and 400 pieces

| Length | colour coded | 100 pcs REF | 400 pcs REF |
|----------|---|----------------|----------------|
| 11 mm |  | 21781 | 21798 |
| 12 mm |  | 21782 | 21795 |
| 13 mm |  | 21783 | 21796 |
| 15 mm |  | 21784 | 21797 |
| Assorted | | 21785 | 21786 |

Articulating paper

Shapes Straight 120 × 20 mm,
horse-shoe 68 × 80 mm,
half horse-shoe* 34 × 80 mm

Thicknesses Thin ~50µ, thick ~80µ

Colours Blue both sides, red both sides,
one side blue and one side red

*As per actual available stock



Straight and Half horse-shoe . Box of 12 books of 12 sheets
Horse-shoe. Box of 6 books of 12 sheets

| Shape | thickness | REF | REF | REF |
|--------------|-----------|-------|-------|-------|
| Straight | thin | 36401 | 36403 | 36405 |
| Straight | thick | 36402 | 36404 | 36406 |
| Horse-shoe | thin | 36411 | 36413 | 36415 |
| Horse-shoe | thick | 36412 | 36414 | 36416 |
| ½ horse-shoe | thin | 36421 | 36423 | 36425 |
| ½ horse-shoe | thick | 36422 | 36424 | 36426 |

Hold

Adhesive for fixing impression materials on impression trays.

Note This product is subject to transport restrictions

| Composition | Spray with gas | | Liquid | Spray without gas |
|-----------------------|----------------|------------|-----------|-------------------|
| | Alginates | Elastomers | Alginates | Alginates |
| Poly-amino-amid Resin | 10 % | — | 15 % | 52 % |
| Silicone Elastomer | — | 18 % | — | — |
| Isopropanol | 25 % | 31 % | 42.5 % | 24 % |
| Toluene | 25 % | 28 % | 42.5 % | 24 % |
| Tetrafluoroethane | 40 % | 23 % | — | — |



| Packaging | | | A REF | B REF |
|-------------------|---------|----------------------------|----------|----------|
| Liquid | bottle | 30 ml | 11130 | — |
| Spray | aerosol | 120 ml (50 ml of adhesive) | 11131 | 11135 |
| Spray without gas | bottle | 100 ml | 11132 | — |

A for Alginates and Hydrocolloids . B for Elastomers

Transparent strips

For composites and other filling materials.
Acetate of cellulose.

Types Straight length 140 mm,
curved length 70 mm
Widths 8 mm, 10 mm
Thicknesses X-thin 0.05 mm, thin 0.075 mm



Envelope of 100 pieces

| Type | thickness | 8 mm REF | 10 mm REF |
|----------|-----------|-------------|--------------|
| Straight | x-thin | 42011 | 42021 |
| Straight | thin | 42012 | 42022 |
| Curved | x-thin | 42211 | 42221 |
| Curved | thin | 42212 | 42222 |

Steel abrasive strips

Long-lasting, sterilisable (dry heat) or disinfected (anti-corrosive solution for rotary instruments), resistant and flexible polishing strips.
Single side (1S) or double side (2S) coated.

Grits Fine, medium, coarse
Length 135 mm
Widths 2, 3, 4, 6, 8 mm
Thicknesses (1S) = 0.15 mm, (2S) = 0.25 mm



Pack of 12 pieces single side 1S (0.15 mm)

| Width | fine REF | medium REF | coarse REF |
|-------|-------------|---------------|---------------|
| 2 mm | — | 46002 | — |
| 3 mm | — | 46012 | — |
| 4 mm | 46021 | 46022 | 46023 |
| 6 mm | 46031 | 46032 | 46033 |
| 8 mm | — | 46042 | — |

Pack of 12 pieces double side 2S (0.25 mm)

| Width | fine REF | medium REF | coarse REF |
|-------|-------------|---------------|---------------|
| 2 mm | — | — | — |
| 3 mm | — | — | — |
| 4 mm | 46121 | 46122 | 46123 |
| 6 mm | — | 46132 | — |
| 8 mm | — | — | — |

Abrasive Polyester finishing strips

End of life

For the finest polishing of Amalgams, Composites, Cements, etc. Aluminium Oxide (Al₂O₃) abrasive, firmly fixed to a resistant and flexible Polyester support. Can be used either wet or dry.

Types White, colour coded
Length 150 mm
Widths 4 mm (medium), 5.5 mm (wide)
Grits X-fine = white, fine = yellow,
medium = red, coarse = blue



Box of 100 strips. Colour coded

| Type | width | x-fine REF | fine REF |
|-------------|---------|---------------|-------------|
| Regular c/c | 4 mm | — | 44422 |
| Regular c/c | 5.5 mm | 44431 | — |
| Regular c/c | assort. | — | — |

Box of 100 strips. White

| Type | width | x-fine REF | fine REF | medium REF | assorted REF |
|---------|---------|---------------|-------------|---------------|-----------------|
| Regular | 4 mm | 44021 | 44022 | — | — |
| Regular | 5.5 mm | 44031 | 44032 | 44033 | 44035 |
| Regular | assort. | — | — | — | — |





FIBRAPOST & SEALACORE

Glass fiber posts with composite FIBRAPOST feature an improved composition (Urethane Dimethacrylate based resin matrix, bis-Phenol A free) with pre-silanisation, higher radiopacity, optimized adhesion to sealing composite and improved mechanical resistance and flexural strength. The SEALACORE family consists in a nano-reinforced self-etch bond and a dual-cure composite cement that are used for cementation of posts and core build-ups.

The SEALACORE product line was specially developed to be used with FIBRAPOST. Thanks to their excellent adhesive properties, SEALACORE products have proved to be exceptionally compatible with FIBRAPOST.

2.2 FIBRAPOST & SEALACORE

| | |
|----------------------|----|
| FIBRAPOST PLUS & LUX | 40 |
| SEALACORE DC | 41 |

FIBRAPOST



FIBRAPOST PLUS

Properties

FIBRAPOST are made from specially treated Glass fibres (approx. 65 %) reinforced with, and firmly bound to, an UDMA (Urethane Dime-thacrylate) resin (approx. 20%) known for its resistance, safety and capacity to bind to Composite resins used for post cementation.

Features and benefits

- Optimal flexural strength (~ 35 GPa) close to the one of the dentine
- High resistance (~ 1'300 MPa) thanks to interlaminar fibres cohesion
- Natural shade for aesthetic reconstruction
- Convenient Radiopacity (~ 200 % Al, Ø 2.5 mm)
- High adhesion capacity between post and cement
- Biocompatibility

Indication

Core build-ups, especially for anterior teeth where aesthetic is of great importance.
To be used with the corresponding DRILL-A-POST, spiral bladed tapered reamers with adjustable depth stopper (page 13)

Packaging

Now also available in hygienic and resistant blisters of 5 posts.
Each post is packaged individually

FIBRAPOST LUX


Combines the features of the FIBRAPOST PLUS with light conduction, both axially (towards the apex) and laterally, by diffusion, thanks to the grooves on the post surface. Used preferably with light cured composites.

Assortment classic

- | | | |
|---|--------------------------------------|---------------------|
| | PLUS REF 22595 | LUX REF 22589 |
| 6 | posts of each N° 1 to 4 = 24 posts | |
| 1 | reamer of each N° 1 to 4 = 4 reamers | |

Assortment in blister

- | | | |
|---|--------------------------------------|---------------------|
| | PLUS REF 22630 | LUX REF 22640 |
| 5 | posts of each N° 1 to 4 = 20 posts | |
| 1 | reamer of each N° 1 to 4 = 4 reamers | |

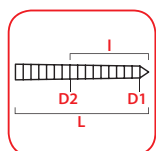

| Refills Sizes of the posts | total length (L) | length of the tape- red part of posts (l) | Ø (mm) at the tip (D1) | Ø (mm) at 11 mm (D2) | colour coded | Packaging classic Rack of 12 pieces | | Packaging blister Box of 5 pieces | |
|----------------------------------|----------------------|--|---------------------------|-------------------------|-----------------|--|------------|--------------------------------------|------------|
| | | | | | | PLUS REF | LUX REF | PLUS REF | LUX REF |
| 1 | 19 mm | 11 mm | 0.85 | 1.30 | □ | 22596 | 22590 | 22631 | 22641 |
| 2 | 19 mm | 11 mm | 1.06 | 1.50 | ■ | 22597 | 22591 | 22632 | 22642 |
| 3 | 19 mm | 11 mm | 1.26 | 1.70 | ■ | 22598 | 22592 | 22633 | 22643 |
| 4 | 19 mm | 11 mm | 1.47 | 1.90 | ■ | 22599 | 22593 | 22634 | 22644 |
| Assort. 1-4 | 3 posts of each size | | | | | 22600 | 22594 | — | — |

SEALACORE DC

Properties

SEALACORE DC is a dual-cure, flowable and adhesive sealing and core build up composite system. It is indicated for the core build up of devitalized tooth as well as for the sealing of glass fiber posts in the root canal. SEALACORE DC is therefore particularly adapted to the seal (or cementation), with one only material, of the FIBRAPOST (PLUS and LUX) in the root canal followed by the core building up

Features and benefits

- Fast hardening by light-curing
- Secure and complete polymerization due to additional chemical curing
- Material savings through optimized mixing syringe
- Low heat development even in voluminous build ups
- Flowable – no plugging
- High adhesion due to self-etch bonding system
- Handy use through application by application tips



SEALACORE DC Bond REF 10471

- | | |
|----|--|
| 1 | dual-cure nano-reinforced self-etch bond (liquid A 4 ml + liquid B 4 ml) |
| 25 | minibrushes |
| 1 | mixing plate |

SEALACORE DC Composite REF 10470

- | | |
|----|------------------------------------|
| 1 | dual-cure composite (syringe 5 ml) |
| 10 | mixers |
| 10 | tips |

| Accessoires | | REF |
|-------------|------------------|-------|
| 50 | minibrushes | 10477 |
| 25 | mixers + 25 tips | 10476 |

ALL-IN-ONE BUILD UP KIT

All you need for a safe and complete core build-up.



ALL-IN-ONE BUILD UP KIT REF 10475

FIBRAPOST PLUS

- | | |
|-------|---|
| 4 × 4 | glass fiber composite posts (sizes 1, 2, 3 & 4) |
| 4 | drills with adjustable depth stopper |

SEALACORE DC

- | | |
|----|--|
| 1 | dual-cure nano-reinforced self-etch bond (liquid A 4 ml + liquid B 4 ml) |
| 1 | dual-cure composite (syringe 5 ml) |
| 10 | mixers |
| 10 | tips |
| 25 | minibrushes |
| 1 | mixing plate |

Goldplated screw posts

Alloy 24 Carat Gold plated Brass

Lengths Short 8 mm, medium 9.5 mm,
long 12 mm, x-long 14 mm

Diameters

N° 1 = Ø 1.05 mm, N° 2 = Ø 1.20 mm


N° 3 = Ø 1.35 mm, N° 4 = Ø 1.50 mm

N° 5 = Ø 1.65 mm, N° 6 = Ø 1.80 mm



Assortment box

| Pieces | keys | REF |
|--------|---------------|-------|
| 60 | 2 | 20701 |
| 120 | 2 | 20752 |
| 240 | 2 | 20700 |
| 1 | cross key | 20756 |
| 1 | socket wrench | 20755 |

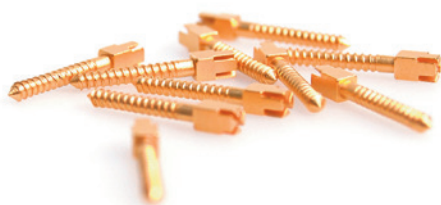

Refill of 12 pieces

| Length | diameter | N°1 Ø1.05 REF | N°2 Ø1.20 REF | N°3 Ø1.35 REF | N°4 Ø1.50 REF | N°5 Ø1.65 REF | N°6 Ø1.80 REF |
|--------|----------|------------------|------------------|------------------|------------------|------------------|------------------|
| Short | 8 mm | 20732 | 20733 | 20734 | 20735 | 20736 | 20737 |
| Medium | 9.5 mm | 20738 | 20739 | 20740 | 20741 | 20742 | 20743 |
| Long | 12 mm | 20744 | 27045 | 20746 | 20747 | 20748 | 20749 |
| X-long | 14 mm | — | — | — | 20750 | — | 20751 |

Refill of 100 pieces

| | | | | | | | |
|--------|--------|-------|-------|-------|-------|-------|-------|
| Short | 8 mm | 20702 | 20703 | 20704 | 20705 | 20706 | 20707 |
| Medium | 9.5 mm | 20708 | 27009 | 20710 | 20711 | 20712 | 20713 |
| Long | 12 mm | 20714 | 27015 | 20716 | 20717 | 20718 | 20719 |

Titanium screw posts

Grooved head for Composites

Lengths and diameters See above



Refill of 12 pieces

| Length | diameter | N°1 Ø1.05 REF | N°2 Ø1.20 REF | N°3 Ø1.35 REF | N°4 Ø1.50 REF | N°5 Ø1.65 REF | N°6 Ø1.80 REF |
|--------|----------|------------------|------------------|------------------|------------------|------------------|------------------|
| Short | 8 mm | 20932 | 20933 | 20934 | 20935 | 20936 | 20937 |
| Medium | 9.5 mm | 20938 | 20939 | 20940 | 20941 | 20942 | 20943 |
| Long | 12 mm | 20944 | 20945 | 20946 | 20947 | 20948 | 20949 |
| X-long | 14 mm | — | — | — | 20950 | — | 20951 |

Assortment Box of 60 pieces + 1 key REF 20901

Special drills and reamers for screw posts

Properties

The sizes (N°) of the instruments correspond to the sizes (N°) of the Gold Plated and Titanium screw posts

RA shank

Lengths 32 mm, 19 mm (working)

Diameters

N° 1 = Ø 0.97 mm, N° 2 = Ø 1.10 mm

N° 3 = Ø 1.25 mm, N° 4 = Ø 1.40 mm

N° 5 = Ø 1.55 mm, N° 6 = Ø 1.70 mm



Drills



Reamers

Rack of 6 pieces

| | N° 1-6 REF |
|---------|---------------|
| Drills | 21007 |
| Reamers | 21017 |



Prophylaxis and Oral Hygiene


| | |
|-------------------|---------------------------------|
| 3.1 | Prophylaxis and Oral Hygiene |
| for Professionals | 46-47 |
| for Patients | 48 |



new

ProphyCups

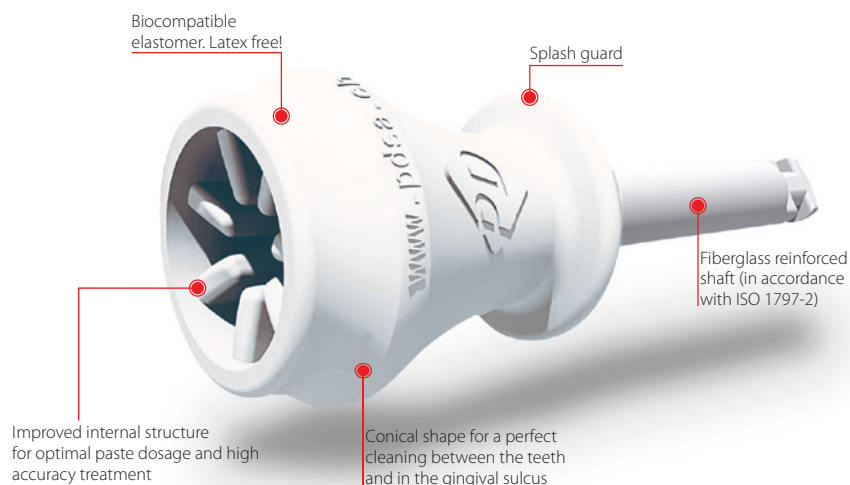Single use prophylactic cups, developed to achieve best cleaning and polishing results by removing all plaques and stains without any risk of scratch marks or damage of the enamel, composite or dentine.



ProphyCup Snow

Presentation

Soft hardness (49 ShA)
for gentle to normal polishing,
cleansing and plaque removal

Packaging

Clean packaged in plastic
pouches of 8 pieces



Packaging

| | |
|------------------|-------|
| Box of 8 pieces | 47051 |
| Box of 48 pieces | 47052 |

REF

ProphyCup Fire

Presentation

Medium hardness (59 ShA)
for medium to heavy
stain and plaque removal

Packaging

Clean packaged in plastic
pouches of 8 pieces



Packaging

| | |
|------------------|-------|
| Box of 8 pieces | 47061 |
| Box of 48 pieces | 47062 |

REF

Rubber polishing cups

Webbed or serrated type for polishing teeth after depuration

Nylon brush type for cleaning and polishing
Mounted on a strong metallic shank, the cups are splash proof and resist autoclaving provided no chemicals are used

Available in 2 grades Soft = red, medium = white



Box of 12 pieces

| Type | soft REF | medium REF |
|---------------------------------------|-------------|---------------|
| RA Webbed | 46502 | 46602 |
| RA Serrated | 46512 | 46612 |
| RA with Nylon brush, Dual-Cups 2 in 1 | 46550 | 46650 |

Prophylactic paste

For cleaning and polishing of all tooth surfaces.

ProphyPaste Pro-N¹⁰⁰

Presentation 100 % natural prophylactic paste

Flavour Mint

Available in 2 grits Fine or medium

Composition Glycerin, calcium carbonate, pumice, deionized water, silica, gelling agent, wetting agent, preservatives, fragrances and dye



| Packaging | fine REF | medium REF |
|---------------|-------------|---------------|
| 3 jars of 50g | 47011 | 47016 |
| 1 jar of 50g | 47010 | 47015 |

ProphyPaste Pro-F^{1.23}

Presentation Prophylactic paste containing 1.23% Fluoride

Flavour Mint

Available in 2 grits Fine or medium

Composition Fluoride, glycerin, calcium carbonate, pumice, deionized water, silica, gelling agent, wetting agent, preservatives, fragrances and dye



| Dispenser box of 125 cups of 2g | REF |
|---------------------------------|-------|
| Fine | 47020 |
| Medium | 47025 |

Pumice powder

Presentation

Quality xxx-fine specially
air screened for dental use

Flavour Peppermint



| Jar | REF |
|-------|-------|
| 500 g | 46403 |
| 1 kg | 46404 |

Effervescent mouth wash tablets

For the preparation of a deodorizing mouth wash and for the effective removal of debris.

Presentation

Pink colour, Peppermint flavoured tablets

Indications

- Chair-side: mouth wash
- At home: daily oral care and hygiene

Composition

1 tablet of 550 mg contains:

Sodium Bicarbonate, Tartaric Acid,
Potassium Bitartrate, Sodium Laurylsulfate,
Magnesium Peroxide, Peppermint Oil, Excipient



| Jar | REF |
|---------------|-------|
| ~ 250 tablets | 66301 |
| ~ 500 tablets | 66302 |
| ~1000 tablets | 66303 |

Disclosing tablets

Indication

Disclosing tablets allow, without
prophylactic or therapeutic action,
to highlight the presence of dental plaque
by means of a red colour food dye

Composition

1 tablet contains:

5 mg (1.36%) Erythrosine (FD&C Red colouring
agent), 2 mg (0.54%) Mint Oil, Excipient ad 360 mg.



| Jar | REF |
|---------------|-------|
| ~ 100 tablets | 11611 |
| ~ 500 tablets | 11613 |
| ~1000 tablets | 11612 |

Dental floss

For the prevention of the dental caries in those interdental spaces which the toothbrush cannot reach.

Types Nylon waxed, pure Silk waxed



Plastic dispenser of 25 m Box of 12 dispensers

| Type | REF |
|-----------|-------|
| Nylon | 66112 |
| Pure Silk | 66012 |

Glass tube with 100 m spool Box of 12 tubes

| Type | REF |
|-----------|-------|
| Nylon | 66113 |
| Pure Silk | 66013 |

Glass tube with 100 m spool

| Type | REF |
|-----------|-------|
| Pure Silk | 66003 |

Dental sticks

Triangular shaped interdental Swiss lime wooden sticks for plaque cleansing and gum stimulation.

Peppermint flavoured for mouth hygiene and clean breath
Handy envelope of 1 plate of 25 snap-off sticks



| Packaging | REF |
|---------------------------------|-------|
| Pack of 12 boxes of 5 envelopes | 65920 |



Special Supplies for the Clinic and the Laboratory

4.1

Special
Supplies

Clinic
Laboratory

52-55
56



Cotton pellets

Soft and consistent in shape, the highly absorbent Cotton pellets are perfect for drying cavities or applying solutions.

5 sizes N°1-5



Box of 6 refills for dispenser model A (see below)

| N° | weight | quantity of pellets | REF |
|----|---------|---------------------|-------|
| 1 | 6 × 3 g | 6 × 105 ≈ 630 | 33251 |
| 2 | 6 × 3 g | 6 × 245 ≈ 1470 | 33252 |
| 3 | 6 × 3 g | 6 × 370 ≈ 2220 | 33253 |
| 4 | 6 × 3 g | 6 × 675 ≈ 4050 | 33254 |
| 5 | 6 × 3 g | 6 × 1035 ≈ 6210 | 33255 |

Box of pellets in bulk for dispenser model B (see below)

| N° | weight | |
|----|--------|-------|
| 1 | 32 g | 33261 |
| 2 | 32 g | 33262 |
| 3 | 18 g | 33263 |
| 4 | 18 g | 33264 |
| 5 | 18 g | 33265 |

Dispenser for Cotton pellets

2 types

Model A for pellets refills

The inside of the dispenser is provided with small clasps which secure the refill cardboard box.

Model B for pellets in bulk

The inside of the dispenser is provided with a spring device which pushes the pellets up.

Each model can be supplied with:

- ▶ Large mesh for pellets N° 1 and N° 2
- ▶ Small mesh for pellets N° 3, N° 4 and N° 5

The dispensers are supplied without pellets.



| Mesh | | model A REF | model B REF |
|-------|-------|----------------|----------------|
| Small | empty | 33210 | 33230 |
| Large | empty | 33211 | 33240 |

EZ Dropper . Dispenser for solutions

The Solution for your Solutions!

A «Must», this elegant, handy multi-dispenser kit for medicaments, disinfectants or chemical solutions should be a permanent feature on the chairtrays in every dental practice.

3 models

Small 3 drop plunger bottles

+ 3 pairs of Dappen dishes

Size of the stand 124 × 80 mm

Regular 4 drop plunger bottles

+ 4 pairs of Dappen dishes

Size of the stand 162 × 86 mm

Improved 3 drop plunger bottles

+ 3 pairs of Dappen dishes

+ 1 Cotton pellet dispenser

model A or B small mesh

Size of the stand 80 × 86 mm

Colour-coded bottles and dishes

Blue, green, yellow-brown, transparent

On request, the bottles are

also available white-opaque

Sterilizable with cold chemicals

Do not autoclave



Complete kit

| Model | mesh | REF |
|------------|-------|-------|
| Small | | 11790 |
| Regular | | 11791 |
| Improved A | small | 11786 |
| Improved A | large | 11787 |
| Improved B | small | 11788 |
| Improved B | large | 11789 |

Drop bottle (15 ml) Per piece

| Colour | transparent REF |
|--------|--------------------|
| Pink | 11760 |
| Yellow | 11761 |
| Blue | 11762 |
| Green | 11763 |

Dappen dish

Glass dish, ideal for holding any type of solution.



Per piece

| Colour | REF |
|-------------|-------|
| Yellow | 11781 |
| Blue | 11782 |
| Green | 11783 |
| Transparent | 11780 |

Bur cleaning brush

Thread Brass wire or Nylon
2 models Metal slide
 or Stainless Steel telescopic handle



Bag of 1 piece

| | REF |
|------------------------------------|-------|
| Metal slide Brass wire | 69302 |
| Metal slide Nylon | 69304 |
| Stainless Steel telescopic handle* | 69303 |

End of life*

Saliva ejectors

Straight, flexible and transparent.

PVC tube with inserted copper plated steel wire
Disposable
Removable PE tips
Lengths 155 mm, 125 mm



Pack of 100 pieces

| Type | length | REF |
|------|--------|-------|
| A | 155 mm | 65311 |
| B | 125 mm | 65312 |

Napkin holders

With 2 «Alligator» clips Nickel plated brass
Autoclavable at 120°C without chemicals



Plastic necklace . Bag of 1 piece

| Colour | length | REF |
|-------------|--------|-------|
| Transparent | 40 cm | 65801 |
| Red | 40 cm | 65802 |
| Blue | 40 cm | 65803 |
| White | 40 cm | 65804 |
| Green | 40 cm | 65805 |
| Yellow | 40 cm | 65806 |
| Violet | 40 cm | 65807 |

Ball chain, chrome plated . Bag of 1 piece

| | | |
|--------|-------------------------------|-------|
| Chrome | regular 40 cm | 65810 |
| Chrome | adjustable length (max 44 cm) | 65811 |

Universal solder, Silver base

Ready to use special Silver base solder, with pre-incorporated flux.

Properties

Can be used to braze many metals/alloys to a like or different metal/alloy as long as their melting points are not lower than ~600°C. Do not use additional flux, cleaning or thinning agents.

This quality is currently used for orthodontics, temporary appliances, repairs, trials, etc.

It offers, in the mouth, the same resistance as 14 Carat Gold Solder.

Melting points 630°C - 680°C

Seamless tube pre-filled with flux

Length 150 mm, Ø 1 mm

White colour (Cadmium free) for Stainless Steel, Chrome Cobalt, Silver, etc.

Yellow colour (contains Cadmium) for Copper, Brass, Gold, etc.



Envelope of 12 sticks

| | REF |
|--------------|-------|
| White | 10001 |
| Yellow* | 10003 |
| End of life* | |

Note End of life*

Soldering tips

End of life

Color Yellow

Without flux

Tip mounted on a Brass shank of 30 mm

Ø 1.10 mm × length 0.85 mm

Ø 1.10 mm × length 1.7 mm



Pack of 100 pieces

| Length | REF |
|---------|-------|
| 0.85 mm | 10060 |
| 1.7 mm | 10061 |

Piezo Gas micro-torch

For orthodontic soldering and brazing.

Temperature of flame 1300 °C

Length of flame adjustable 10 to 35 mm

Time span for use 25 min



Complete device

REF 10036

This image shows a single sheet of white paper with horizontal ruling lines. The lines are evenly spaced and run across the width of the page. There are no margins, text, or other markings on the paper.

Handwriting practice lines consisting of 30 horizontal rows. Each row is defined by two thin grey lines, with a slightly thicker grey line in the middle of each row to guide letter height.

A-B

| | |
|--|----|
| Abrasive Polyester finishing strips..... | 36 |
| Acrylic points..... | 8 |
| Aluminium crowns, preformed..... | 32 |
| Articulating paper..... | 35 |
| Bur cleaning brush..... | 54 |

C

| | |
|--|----|
| « Cavity Liner » | |
| Calcium Hydroxide paste..... | 20 |
| Calcium Hydroxide paste, light-curing..... | 20 |
| Calcium Hydroxide compound..... | 21 |
| Calcium Hydroxide powder..... | 21 |
| C&B Zinc Phosphate cement..... | 24 |
| Copal varnish..... | 22 |
| Copper bands..... | 30 |
| Cotton pellets..... | 52 |

D

| | |
|---|----|
| Dappen dish..... | 53 |
| Dental floss..... | 49 |
| Dental sticks..... | 49 |
| Disclosing tablets..... | 48 |
| Dispenser for cotton pellets..... | 52 |
| DRILL-A-POST . Spiral bladed tapered reamers..... | 11 |

E

| | |
|---|----|
| EDTA 17%..... | 17 |
| « Endo » plastic pipettes..... | 16 |
| Eugenate Desobturator..... | 17 |
| Eugenol..... | 27 |
| EZ Dropper . Dispenser for solutions..... | 53 |

F**FIBRAPOST**

| | |
|--|----|
| FIBRAPOST PLUS & LUX..... | 40 |
| FIBRAPOST & SEALACORE . All-in-one Build Up Kit..... | 41 |

G-H-J

| | |
|---|----|
| Gates-Glidden drills..... | 10 |
| Glass Ionomer cements | |
| Ionover Filling . Glass Ionomer Filling cement..... | 26 |
| Ionover Luting . Glass Ionomer Luting cement..... | 26 |
| Gutta-percha points..... | 8 |
| Gutta-percha stopping..... | 23 |
| Hold..... | 35 |
| Jiffy tubes for cement..... | 27 |

M-N**MAP SYSTEM**

| | |
|-------------------------|----|
| Intro kit NiTi..... | 14 |
| Universal kit NiTi..... | 14 |
| Surgical kit..... | 14 |
| PD MTA White..... | 15 |
| Needles..... | 15 |
| Accessories..... | 15 |

| | |
|---------------------------------------|----|
| Marking rings..... | 16 |
| Mouth wash tablets, effervescent..... | 48 |
| Napkin holders..... | 55 |

P

| | |
|---|----|
| Paper points..... | 9 |
| Peeso reamers..... | 10 |
| PELLA transparent crown forms..... | 28 |
| Piezo Gas micro-torch..... | 56 |
| Plastic pipettes | |
| « Endo »..... | 16 |
| « One-piece »..... | 16 |
| Polishing cups, rubber..... | 46 |
| Prophylactic paste..... | 46 |
| ProphyPaste Pro-F ^{1,23} | 46 |
| ProphyPaste Pro-N ¹⁰⁰ | 46 |
| Pro-Pul-Pan . Zinc Oxide-Eugenol reinforced cement..... | 24 |
| Pumice powder..... | 48 |

R

| | |
|--|----|
| Resin Remover..... | 17 |
| Root canal barbed broaches..... | 11 |
| Rosen root canal cement . Rickert's paste..... | 25 |
| Rosen root canal cement Silver free..... | 25 |

S

| | |
|---------------------------------|----|
| Saliva ejectors..... | 54 |
| Screw Posts | |
| Goldplated..... | 42 |
| Titanium..... | 42 |
| Special drills and reamers..... | 43 |

SEALACORE

| | |
|--|----|
| SEALACORE DC..... | 41 |
| SEALACORE DC Bond..... | 41 |
| SEALACORE DC Composite..... | 41 |
| SEALACORE & FIBRAPOST . All-in-one Build Up Kit..... | 41 |
| Soldering tips..... | 56 |
| Steel abrasive strips..... | 37 |
| Stripper, adjustable..... | 37 |

T-U-V-W-Z**TEMPOTEC**

| | |
|---|----|
| Temporary Eugenol free filling material in paste..... | 23 |
| Tofflemire matrix bands..... | 33 |
| Tofflemire matrix retainers..... | 33 |
| Transparent rings..... | 31 |
| Transparent strips..... | 36 |
| Universal solder, Silver base..... | 56 |
| Vaseline..... | 22 |
| Wooden Wedges | |
| Interdental..... | 34 |
| Anatomical..... | 34 |
| Zinc Oxide..... | 27 |
| Zinc Oxide-Eugenol reinforced cement . Pro-Pul-Pan..... | 24 |



Produits Dentaires SA
Vevey . Switzerland

Produits Dentaires SA
Rue des Bosquets 18
1800 Vevey
Switzerland
Phone +41 21 921 26 31
Fax +41 21 921 39 79
info@pdsa.ch
www.pdsa.ch